

Basic Details & Measurements

Scan

Locking Mechanism

Modifications

Finishing

Basic Details

Patient Name *

Patient Name

Date of Birth *

dd/mm/yyyy

Height (cm)

65

Weight (kgs) *

50

Mobile Number

10 digit phone number

Email

Email

Gender *

Select

Amputation Date

dd/mm/yyyy

Amputation Leg

Select

Reason of Amputation

Select

Activity Level *

Select

Socket Type *

Check Socket

Design Variation *

Standard (sx)

Model *

AddEase

Measurements

A

C

Patela Tendon

Patela Tendon to Ground (cm)

A: Stump Length (cm)

B

Ground

A Stump Length

Sound Leg Length (Patela Tendon to Ground)

Circumference of Stump at 5 cm gap

Short Stump

Upto 18 cm

Medium Stump

19-24 cm

Long Stump

25-32 cm

Very Long Stump

Above 32 cm

Value A Stump Length (cm) *

20

Value B Stump Size (cm)

20

Foot Type

Select

Shoe Size (cm)

65

Flexion Angle(Deg)

(Deg)

Add/Abd Angle (Deg)

(Deg)

Stump Type

Select

Left Anterior


Left Lateral

Left Posterior

Stump Condition (please describe any specific condition of the stump example bony prominence etc.)

Previous Prosthetic Experience (Please describe any previous experience of Prosthetics used, Make, Model, Type, Issues with it and expectation from the new Prosthetic socket)

Change image and add table for C

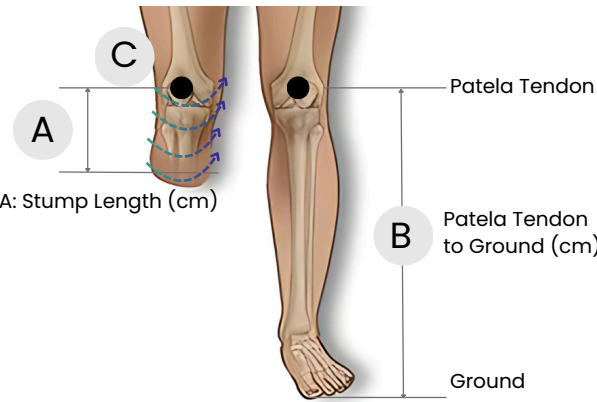

- A Stump Length
- B Sound Leg Length (Patela Tendon to Ground)
- C Circumference of Stump at 5 cm gap

Short Stump	Upto 18 cm	Medium Stump	19-24 cm
Long Stump	25-32 cm	Very Long Stump	Above 32 cm

Make of equal length

Value A Stump Length (cm) *

20

Value B Stump Size (cm)

20

Value C - Circumference of Stmp at 5 cm gap (cm)

0 cm		15 cm	
0 cm		20 cm	
5 cm		25 cm	
10 cm		30 cm	

Socket Selection Criteria

Standard (SX)

Socket with Core
functionality and durability



Adaptive (AX)

Adjustable Socket for varying
residual limb conditions



AddiEase

Premium Definitive Sockets
Printed on HP-MJF



AddiEaseEco

Economy Definitive & Check Sockets
Printed on FDM Printer



AddiEaseMould


Moulds
Printed on FDM Printer



AddiEaseMould-HR

Heat Resistant Moulds
Printed on FDM Printer



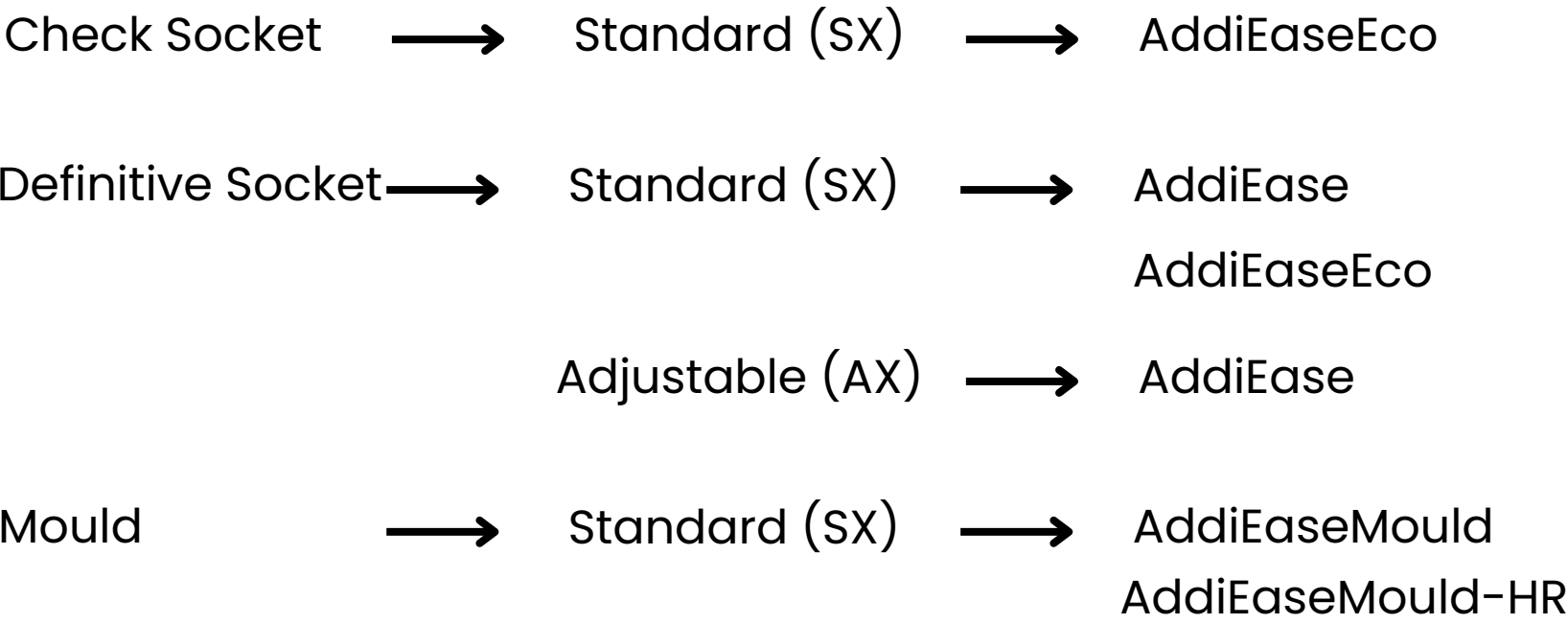

Scan

Scan Condition

Direct Body

With Liner

Liner Thickness

1 mm 2 mm 3 mm 6 mm 8 mm

Liner Tpe

Socks Cushion Locking Seal-In

Scan Upload

Upload Scan

Both ▾

Left Foot

Right Foot

Upload Additional Files
(Design / Rough calculations etc.)

File 1

File 2

Upload Link with Photos
Upload in Google /Cloud drive and give relevant permission

[https://drive.google.com/...](https://drive.google.com/)

Locking Mechaism

Locking/Suspension System

- Bull Dog
- Manual Suction
- Auto Expulsion Suction
- Fit Kit Lock/ Shuttle Lock
- Seal-in Suction

- None
- Other
- Provision for Dummy

Provide Clickable Links to open a pdf file attached with each

Note: 4 hole adapter shall be provided


Add File with dimensions

Design Modications

Delete
this

Socket Design

Please specify the design considerations for each point from A to N. Use "-" to indicate Apply pressure (Reduction) and "+" to indicate Relief at the particular area. All values should be in millimetres (mm). For eg for applying reduction of 6 mm at Patela Tendon, please specify -6

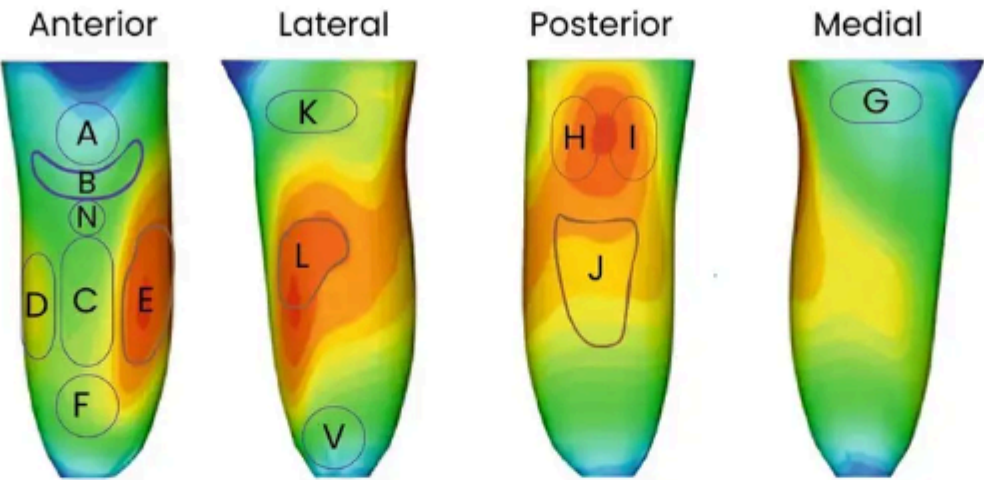

Markings on the Residual Limb prior to Scan & Pressure/Relief areas

Global Volume Reduction
(please specify the percentage reduction in Volume without reducing the length of the socket)

Default Value:2%

Validation – Min 0, max :5%

Please note as a general guideline our design algorithm will consider an overall 2% reduction in the stump dimensions to design the Socket. The below values for "+" and "-" should be done based on this assumption.

Delete
this

A. Patella (mm) Default Value +1	B. Patella Tendon (mm) Default Value -8	C. Crest of tibia (mm) Default Value +3
D. Lateral Shaft Tibia (mm) Default Value -3	E. Medial Shaft Tibia (mm) Default Value -3	F. Distal End Tibia (mm) Default Value +6
G. Medial Supracondylar (mm) Default Value -4	H. Hamstring Tendons (mm) Default Value +2	I. Hamstring Tendons (mm) Default Value +2
J. Posterior Connpartment (mm) Default Value -4	K. Lateral Supracondylar (mm) Default Value -4	L. Head of Fibula (mm) Default Value +3
M. Distal End Fibula (mm) Default Value +4	N. Tibial Tuberosity (mm) Default Value +1	

WhatsApp

Submit

Finishing

Design & Printing

Design by :

Self Addiwise

Print by :

Self Addiwise

Extras For AddiEase

Latices

Yes No

Finish :

For AddiEase



Bead Blasting



Black Dye



Colour

For AddiEaseEco



Bead Blasting



Black Dye



Epoxy

For AddiEaseMould/AddiEase Mould-HR

Not Applicable

If Design by Self, and Print by Self

Sales Order is Created at Zero Cost

If Design by Self, and Print by Addiwise

AddiCoin shall be consumed for Design and Sales Order is Created – Print order

If Design by Addiwise and Print by Addiwise

Two Suborders will be created – One for Design and One for Print

Integration with
Payment Gateway

- Design Sub Order will pick price from Design Price List
- Print Sub Order will pick price from Print Price List



Dashboard

Orders

Organization

Products

Prosthetics Lower Extremity



Below Knee (BK)



Above Knee (AK)

Prosthetics Upper Extremity



Below Elbow (BE)



Above Elbow (AE)



Shoulder Disarticulation



AddiStud

Orthotics



Insoles



AFO/SMO



Braces



Cranial Helmet



Face Mask

Scanners



AddiPrime



AddiPrime RP

Printers



HP-MJF



AddiPrint 1



AddiPrint 2

AddiNxT Design Platform



Standard



Premium